

# Compat ella<sup>®</sup>

## Compat Ella<sup>®</sup> Dual Spike Set

The Compat Ella<sup>®</sup> Dual Spike Set is intended for the simultaneous administration of enteral nutrition and/or hydration via two separately spiked containers and is suitable for use with the Compat Ella<sup>®</sup> feeding pump.

### Features at a Glance

The Compat Ella<sup>®</sup> Dual Spike Set features two ENPlus nutrition connectors allowing **simultaneous enteral feeding and/or hydration**, an ENFit<sup>®</sup> compatible 3-way medication port and an ENFit<sup>®</sup> distal end (patient side tube connection).


**Reliable Feeding:** The Compat Ella<sup>®</sup> Dual Spike Set features patented fail-safe Auto-Anti-Free-Flow (AAFF) technology. AAFF is incorporated into the cassette and is always automatically activated to **ensure patient safety**.

**Safer Connectors:** The Compat Ella<sup>®</sup> Dual Spike includes ENFit<sup>®</sup> connectors made according to ISO 80369-3 and two ENPlus connectors made according to ISO 18250-3. ENFit<sup>®</sup> and ENPlus connectors are designed to be IV incompatible to **promote patient safety** by helping to avoid accidental misconnections.

**Safe Materials:** The Compat Ella<sup>®</sup> Dual Spike Set is sterile and **DEHP, BPA and Latex-Free\***.

#### Compatibility:

- Nutrition containers: ENPlus
- Enteral feeding tubes: ENFit<sup>®</sup>
- Syringes: ENFit<sup>®</sup>

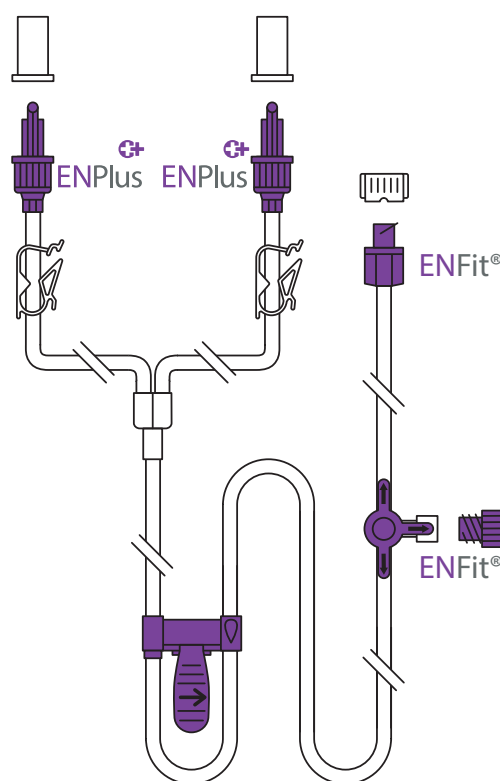


### Ordering Information

Product Name	Product Code	Qty/Box
Compat Ella <sup>®</sup> Dual Spike Set	12297288	30

*Have you met Compat Ella<sup>®</sup>?*

# Compat Ella® Dual Spike Set

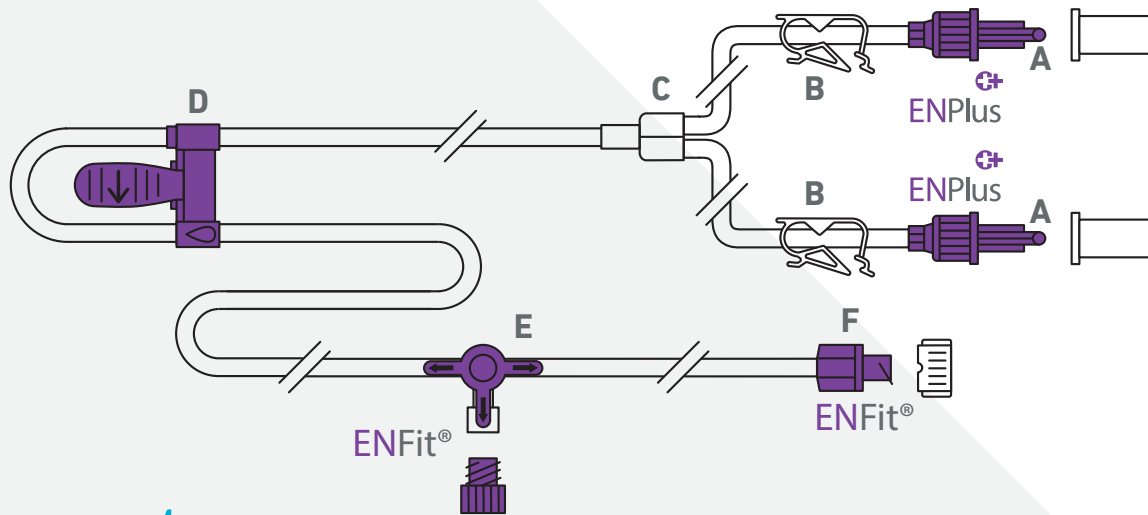

## Components

- A ENPlus 'enteral only' connectors (nutrition connection) with revolving lock and cap
- B Pinch clamps to manually shut off the flow of fluids
- C Y junction
- D Compat Ella® pump cassette with integrated and patented Auto-Anti-Free-Flow (AAFF) safety feature
- E ENFit® compatible 3-way medication port with cap
- F ENFit® distal end with cap (patient side tube connection)

## Specifications

- Sterile:** Sterilized with Ethylene Oxide (EO).
- Packing:** Individually packed in cases of 30 pieces (29.5 x 26 x 19 cm).
- Storage:** Store at room temperature. Keep away from sunlight and humidity.
- Shelf life:** 3 years from production date. Expiry date on the package.
- Duration of use:** 24 hours. Single Use.
- Length:** Distance from nutrition connection to patient side tube connection: 230 cm.  
Distance from nutrition connection to Y junction: 30 cm.  
Distance from nutrition connector to pump: 50 cm.
- Diameters:** Nominal inner diameter: 3.0 mm. Nominal outer diameter: 4.15 mm.
- ISO Standard:** ENFit® connector made according to ISO 80369-3 and ENPlus connector made according to ISO 18250-3.
- CE mark:** Complies with the European Medical Devices Directive 93/42/EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
- Classification:** Class IIa.
- Legal manufacturer:** Cedic Srl, Via Liberazione 63/9, IT-20068 Peschiera Borromeo (MI).
- Sustainability:** Environmental Management System complies with ISO 14001:2015.
- Materials:** Tubing: PVC and silicone. All materials DEHP, BPA and Latex-free\*. Additional details on request.
- Country of origin:** Made in Italy.

\*DEHP and BPA are not intentionally added during the manufacturing process. Not formulated with natural rubber latex.

To obtain further information please contact your Compat® representative or go to [www.compat.com](http://www.compat.com)

®Unless stated otherwise, all trademarks are used with permission.

AXI00010 REV.1 16 MARCH 2021

# Compat®